



Graver Technologies

FILTRATION | SEPARATION | PURIFICATION



ZTEC™ P Series Filter Cartridges

*Pleated Polyethersulfone (PES)
Membrane for Sterile Filtration*

Product Specifications

Media: Asymmetric
Polyethersulfone Membrane

Inner core, end caps, cage: Polypropylene

Support layers: Spunbonded
Polypropylene

O-Rings:

Buna-N, EPDM, Silicone, Teflon
Encapsulated Viton O-Rings, Viton

O-Ring Insert: PBT

Micron rating: 0.2 µm

Dimensions

Nominal lengths:

10" 20" 30" 40"

25.4 50.8 76.2 101.6cm

Outside diameter: 2.7" (6.9 cm)

Inside diameter: 1.0" (2.54 cm)

Surface area: 6.8 ft² (0.63 m²)
per 10" element


Operating Parameters

**Maximum sustained
operating temperature:**
176°F (80°C) at 20 psid (1.38 bar)

Maximum differential pressure:
80 psid @ 70°F (5.5 bar @ 21°C)
40 psid @ 160°F (2.8 bar @ 71°C)

Maximum reverse differential pressure:
40 psid @ 70°F (2.8 bar @ 21°C)

Recommended change-out pressure:
35 psid (2.4 bar)

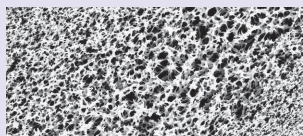

ZTEC P Sterilizing Grade membrane cartridges are validated for complete bacterial retention to yield product sterility in biopharmaceutical final filtration applications. The naturally hydrophilic and low protein binding characteristics of polyethersulfone membrane ensure maximum transmission of active ingredients making it ideal for a wide range of pharmaceutical and biological liquid applications, including the filtration of therapeutics, vaccines, antibiotics, bulk pharmaceutical and other critical biotechnology products. The double-layer PES 0.2 micron membrane filters are manufactured in a cleanroom environment, and integrity tested before shipment to assure consistent performance and quality.

FEATURES & BENEFITS

- Manufactured in an ISO Class 7 Cleanroom Environment
- 100% flushed with ultrapure DI water
- Meets ASTM Standards for Sterility
- Repeatably Steamable/Sanitizable
- 100% Integrity tested prior to release
- Pore size, lot and serial number are stamped on each filter element for identification and traceability
- Complete validation guide available
- Every cartridge comes standard with an embedded O-ring support ring

CERTIFICATIONS

- USP Class VI: Meets USP Class VI Biological Test for Plastics
- FDA Listed Materials: All materials comply with FDA Title 21 of the Code of Federal Regulations Sections 174.5, 177.1520 (polypropylene), and 177.2440 (PES Resin) as applicable for food and beverage contact.
- European Directive for Food Contact: European Regulation No. 1935/2004 and European Regulation 10/2011: Tested for migration behavior and is suitable for contact with all kinds of foodstuffs with minimal rinse-up. Data available upon request.

TYPICAL APPLICATIONS

- Diagnostics
- Ophthalmic Solutions
- Culture Media
- Reagent Chemicals
- LVPs
- Buffers
- Vaccines

PERFORMANCE SPECIFICATIONS


- Hot DI Water: Filter cartridge will withstand temperatures of 185°F (85°C) for up to 30 consecutive minutes.
- Cleaning/Sanitization: Compatible with most common chemical cleaning, sanitizing and sterilizing agents and with pH range from 1-14. Consult factory for specific compatibility information.
- Steam/Autoclave: Cartridges may be steamed or autoclaved for at least 50 thirty-minute cycles @ 275°F (135°C).
- Typical Bacterial Retention Performance: Cartridges have been validated for the complete retention of *Brevundimonas diminuta* at a challenge level of 10⁷ organisms/cm² as prescribed in ASTM 838-05.

ZTEC P NOMENCLATURE INFORMATION

| Filter Type | Retention Rating (microns) | Nominal Length (inches) | End Configuration | Gasket or O-Ring |
|----------------------------|----------------------------|------------------------------------|--|--|
| ZTEC P Series | 0.2 | -5 -30 -10 -40 -20 | P2 226/Flat Single Open End P3 222/Flat Single Open End P7 226/Fin Single Open End P8 222/Fin Single Open End | B Buna-N E EPDM S Silicone T Teflon encap. Viton (O-Rings only) V Viton |
| Example: ZTEC P 0.2-20 P2S | | | | |
| ZTEC P | 0.2 | -20 | P2 | S |

ZTEC P FLOW RATE

Typical Flow Rate Clean Water at Ambient Temperature (per 10" cartridge)

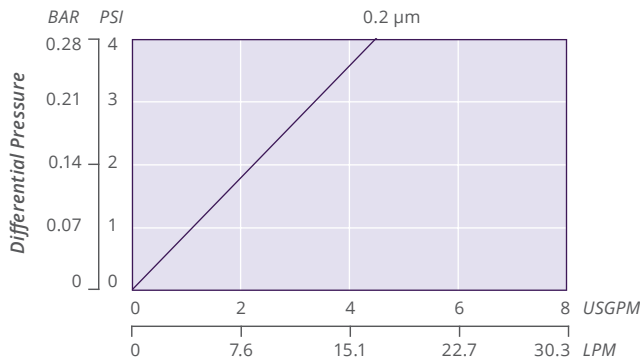

For liquids other than water, multiply pressure drop by the fluid viscosity in centipoise

INTEGRITY TEST SPECIFICATIONS

Maximum Diffusive Air Flow (per 10-inch cartridge) values for ZTEC P filters wet with water:

| Pore Size | Bubble Point | Diffusive Air Flow |
|-----------|---------------------|---------------------------------|
| 0.2 µm | ≥ 40 psig (2.8 bar) | ≤ 30 cc/min @ 32 psig (2.2 bar) |

FOR MORE INFORMATION

Customer Service/Technical Support: 1-888-353-0303
China: +86-21-5238-6576 Asia: +65-9671-9966

GTX-334 3-24

DISTRIBUTED BY

All information and recommendations appearing in this bulletin concerning the use of products described herein are based on tests believed to be reliable. However, it is the user's responsibility to determine the suitability for his own use of such products. Since the actual use by others is beyond our control, no guarantee, expressed or implied, is made by Graver Technologies as to the effects of such use or the results to be obtained. Graver Technologies assumes no liability arising out of the use by others of such products. Nor is the information herein to be construed as absolutely complete, since additional information may be necessary or desirable when particular or exceptional conditions or circumstances exist or because of applicable laws or government regulations. TPE is a trademark of Graver Technologies, LLC.



Graver Technologies | 302-731-1700 | 800-249-1990
info@gravertech.com | www.gravertech.com

A member of The Marmon Group—A Berkshire Hathaway Company